


HOWARD WRIGHT 
PREMA STRETCHER

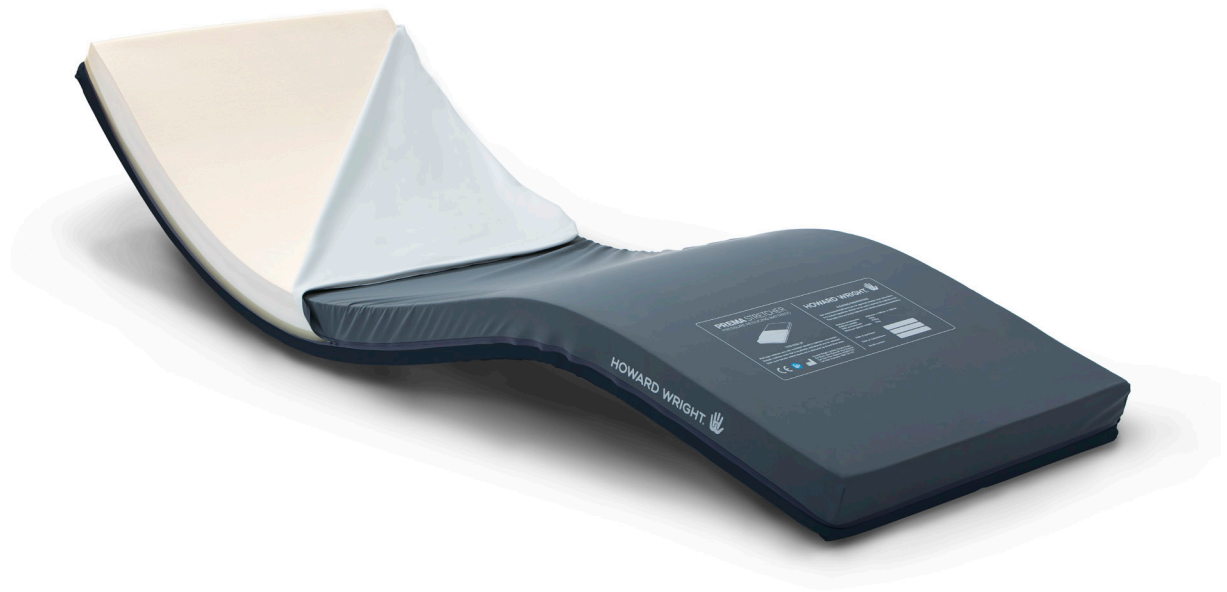

A superior cell structure for uniform head to toe support.

The PREMA Stretcher (Pressure REducing MAttress) and PREMA Stretcher 125 offers patients enhanced levels of comfort and care from its superior cell structure. The mattress has a rotatable design to maximise its service life.


The waterproof and bielastic outer cover is breathable, antimicrobial and fire retardent. All corner seams are welded and a waterproof zip allows easy maintenance.

The materials used are resistant to a range of cleaners and disinfectants including 0.1% sodium hypochlorite. Asset management and auditing are made easier using the on-product information annotation areas.

FEATURES

- High quality foams with superior cell structure
- High flame resistance
- Waterproof, pliable, breathable and antimicrobial outer cover
- Welded seams with a waterproof zip around 3 sides
- X-ray radiolucent foam and cover materials
- Non-slip mattress base material
- Mattress top cover tested to Oeko-Tex Standard® 100
- Mattress top cover halogen free
- 100mm or 125mm thickness available



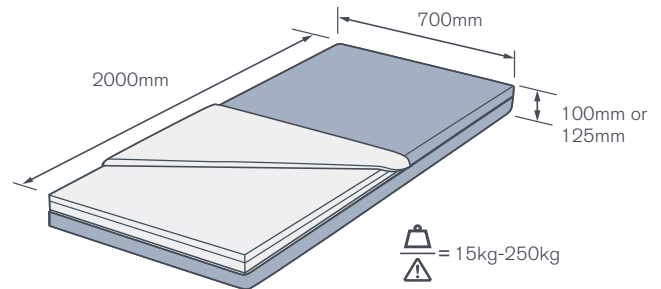

SPECIFICATIONS

MODELS:	PREMA STRETCHER PREMA STRETCHER 125
PART NUMBER:	HCS100C (non slip mattress base) HCS105 (non slip mattress base)
MAXIMUM SWL:	250kg
MATTRESS MASS:	10kg (HCS100C) 13kg (HCS105)

COVER MATERIAL

TOP MATERIAL	BASE MATERIAL
• Material type: bielastic	bielastic
• Composition: 62% PA, 38% PU	62% PA, 38% PU by weight
• Stretch: 2 way	2 way
• Tear strength: 55N Warp, 60N weft	90N warp, 90N weft
• Removable: yes	yes
• Hydrostatic head test: >2M	>2M class DIN53.866
• Water resistant: yes	yes
• Bacteria and mite resistant: yes	yes
• Urine and blood resistant: yes	yes
• Air and vapour permeable: yes	yes
• Cleaning: yes	yes

To wash:

- Wipe with a sponge or soft cloth wetted with warm water containing soap or mild detergent
- Cover can be machine washed in 95°

To disinfect:

- Wash the surface with one of the approved cleaners
- 0.1% sodium hypochlorite solution can be used to disinfect the cover
- Wipe off residue with water using a soft cloth or sponge
- Thoroughly dry
- See our website for detailed information of approved cleaners

FOAM CORE

- Number of layers: 2
- Surface profile: flat
- Density: 50/38
- Hardness: 120/200
- Cleaning: not to be cleaned

MATTRESS COMPLETE

- Skin friction & shear risk: low
- Fire retardant Crib 5 - BS7175
- Pressure mapping for HCS100C. <30.5mmHg. Average 26.25mmHg over 62kg - 128kg weight range
- Pressure mapping for HCS105 - <30.1mmHg. Average 25mmHg over 62kg - 128kg weight range
- Comfort: high
- Anti static: no
- Directional: no
- Sides for use: 1

DIMENSIONS

- PREMA Stretcher Mattress (HCS100C) : L:2000mm, W:700mm, D:100mm
- PREMA Stretcher 125 (HCS105): L:2000mm, W:700mm, D:125mm

www.howardwrightcares.co.uk



Design and Development, Production, Sales, Distribution and Servicing of Medical Beds, Stretchers and Accessories

Howard Wright Limited has a policy of continuous product improvement and reserves the right to change product specifications without notice.
Howard Wright Limited is certified to ISO13485:2012